
CONTACT ARTIKULATOR

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Instructions for use

Mode d'emploi

Istruzioni d'uso

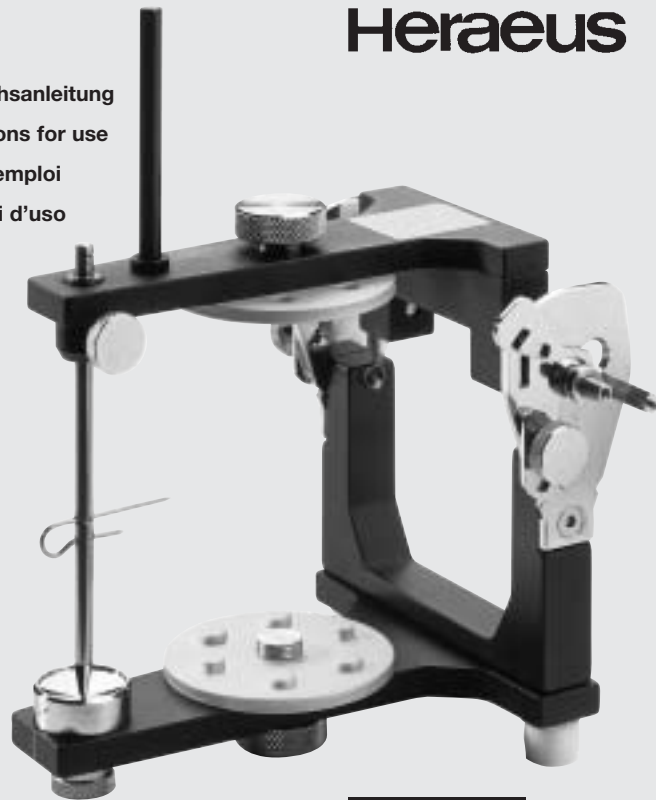Heraeus

**Heraeus
Kulzer**

Partnership-First.

Das Contact System

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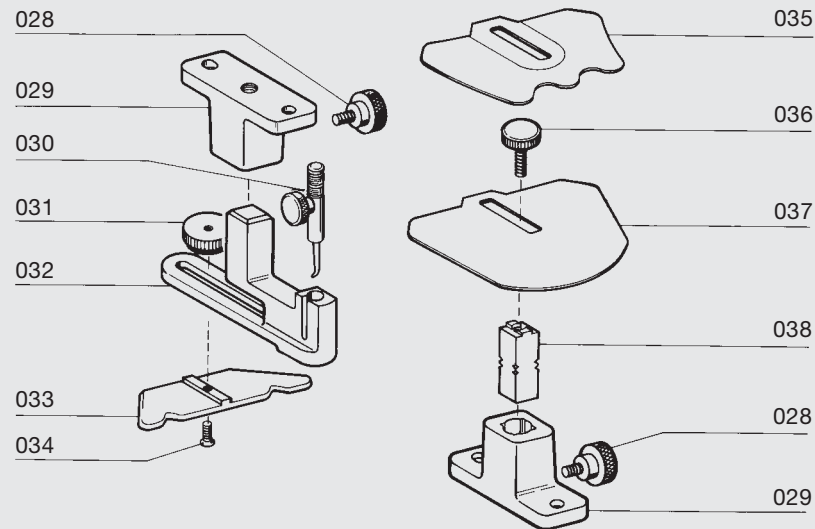
The Contact System

■ The production of prosthetic work that will function correctly and satisfy the patient requires the use of an articulator that carries out the same physiological movements as the patient. For the majority of prosthetic work a semi-adjustable articulator with a number of additional features is completely adequate.

■ The Heraeus Kulzer articulator meets these demands. The Contact System includes registration blocks for intraoral pin

registration methods. These are positioned according to the anatomical and physiological factors and offer the best possible conditions for a reproducible registration and transfer of the corresponding functions of the natural mastication.

■ A wide range of accessories is available for the Heraeus Kulzer Contact Articulator System.

■ To fit all arbitrary face bows, the arbitrary pins supplied are inserted into the corresponding openings in the articulator joints.

■ A split cast system is available as an accessory.

Contact - An Articulator with special quality characteristics.

Semi-adjustable articulator

Our Contact semi-adjustable articulator provides the following adjustments:

Condylar path of 30° relative to the occlusal plane. The distance from the articulator joints to the incisal pin is 105 mm, according to the Bonwill triangle. The inclination of the incisal plate is 15°. If required, an incisal plate with 0° is available, as an accessory. The space available for mounting models is 100 mm.

The condylar joints have a double cone shape and correspond to Prof. Gerber's requirements. The inside cone, with an angle of 17°, and the outside cone, with an angle of 12°, execute the Bennett movement.

Position A

- Stop plate down -

For mounting models and for all opening and closing movements from centric (check of occlusion).

The zero mark on the incisal pin is used to indicate the central occlusal plane.

Position B

- Stop plate up -


For the joint-like reproduction of protrusive and laterotrusive movements.

Position C


- Stop plate centre forward -

For checking retrusion.


The following movements can be executed by means of the joint locking device:

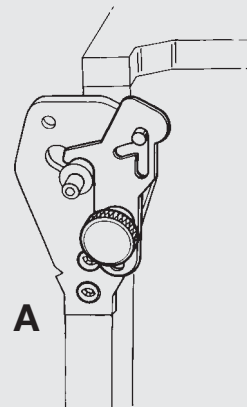A

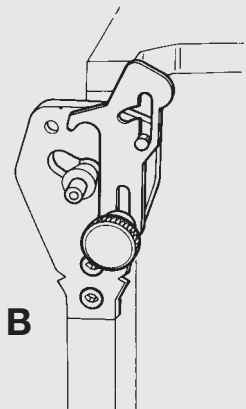
B

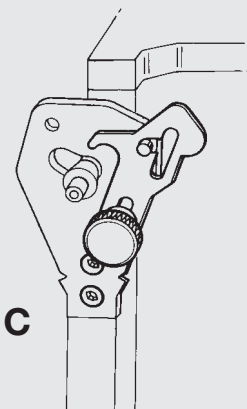C

Accessories

- 1 Adapter for:
- setup calottes
- mounting table
- setup key
- 2 Setup calotte
- 3 Mounting table
- 4 Setup key with sliding plate and symphysis indicator
- 5 Pin support
- 6 Spacer, short
- 7 Spacer, long
- 8 Mounting plates (2 pcs)
- 9 Incisal indicator
- 10 Incisal pin, straight
- 11 Incisal plate 0°, 15°
- 12 Base plates (2 pcs)
- 13 Spacer set, for use with long straight, incisal pins.
- 14 Split cast system on demand

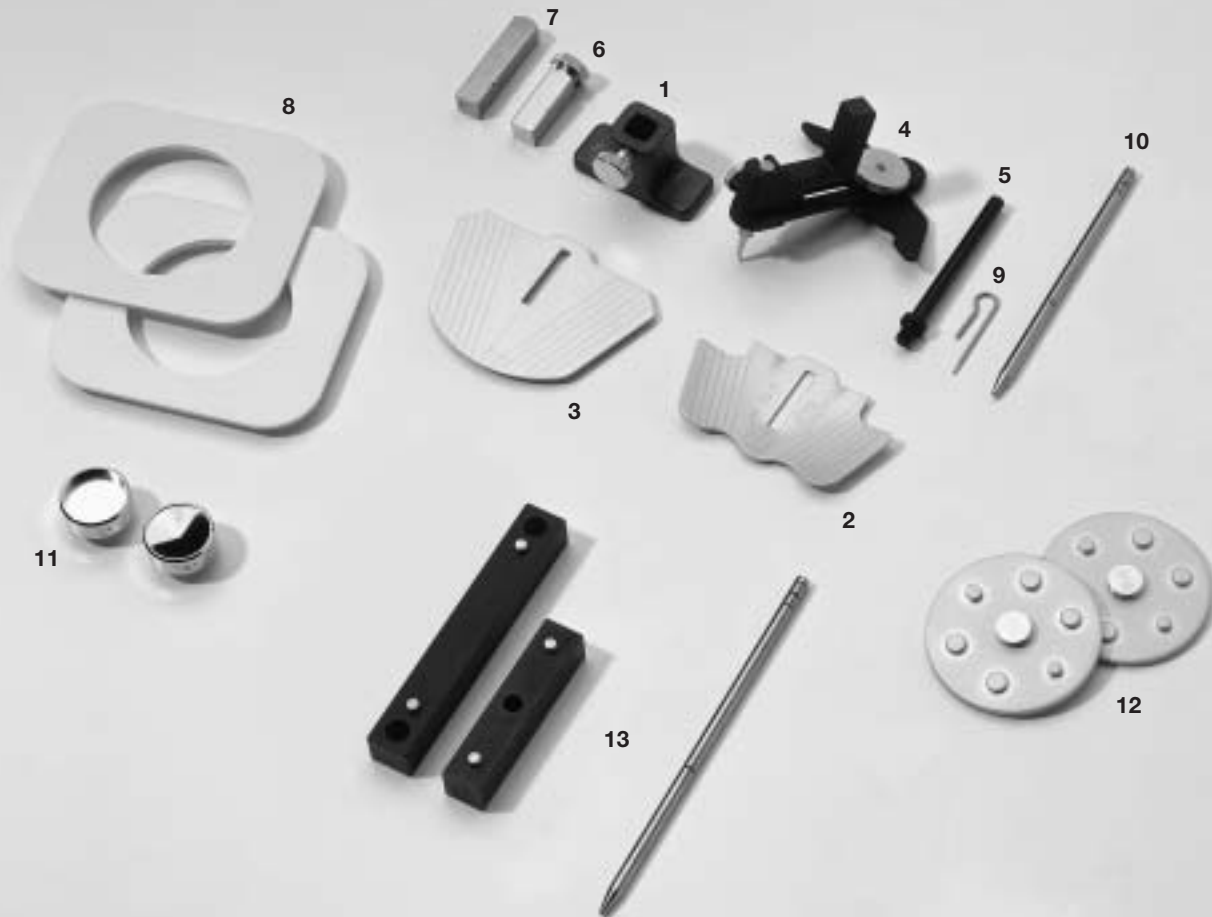
Spare Parts List

- 01 Top section
- 02 Upper mounting
- 03 Joint
- 04 Left joint casing with positioning holes for arbitrary pins
- 05 Positioning screw for joint
- 06 Set-screw for upper mounting
- 07 Lower mounting
- 08 Mounting plate
- 09 Plaster mounting aid
- 010 Bottom section
- 011 Washer
- 012 Set-screw for lower mounting
- 013 Cap
- 014 Positioning screw for mounting plate
- 015 Support pin for top section
- 016 Positioning screw for incisal pin
- 017 Right joint casing with positioning holes for arbitrary pins
- 018 Set-screws for joint casing


- 019 Tension spring
- 020 Locking device for joint
- 021 Washer
- 022 Positioning screw for joint locking device
- 023 Incisal pin
- 024 Incisal positioning pin
- 025 Securing ring
- 026 Incisal plate
- 027 Positioning screw for incisal plate

In case of defect parts, or missing parts please order according to the diagram and the spare parts list, and state the number of the part.

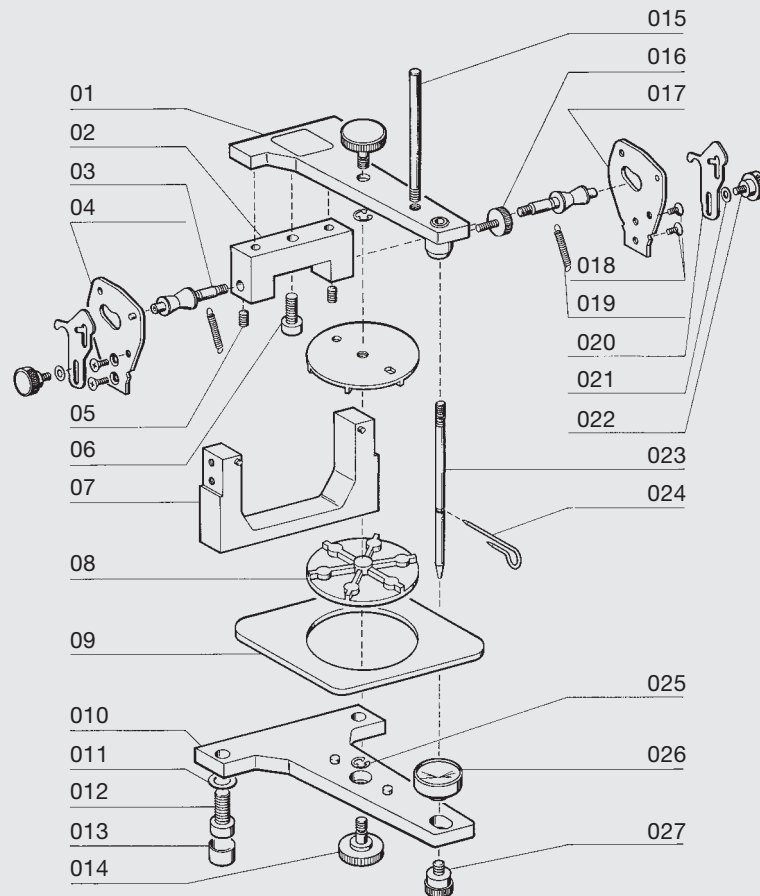Accessories- Spare Parts List

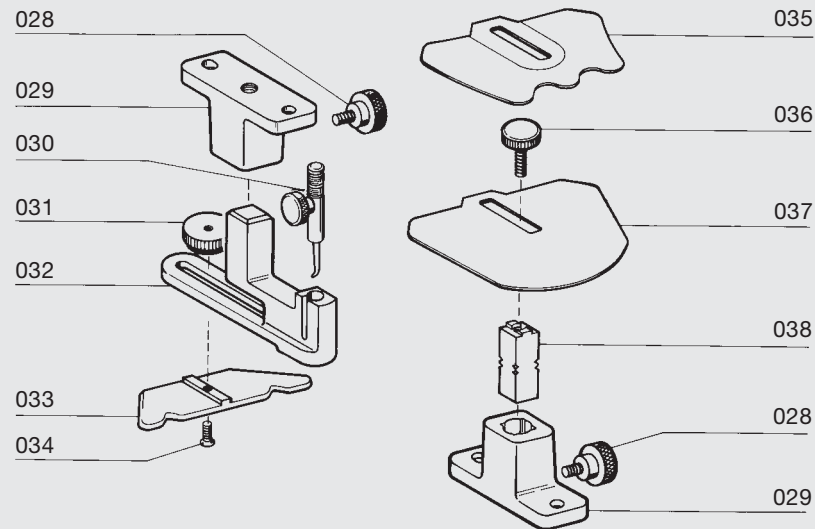- 028 Positioning screw for adapter
- 029 Adapter
- 030 Symphysis indicator with positioning screw
- 031 Positioning screw for support plate
- 032 Mounting key
- 033 Plate for mounting key
- 034 Counter nut for support plate
- 035 Setup calotte
- 036 Positioning screw for calotte and mounting table
- 037 Mounting table
- 038 Distance guage for adapter

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Heraeus Kulzer GmbH & Co. KG